

# Spiegelberg:

Technology for brains




IAP-Catheter for  
Intraabdominal  
Pressure Monitoring

# IAP-Monitoring

The IAP-Catheter monitors Intraabdominal Pressure continuously. It is inserted like a naso-gastric tube. It is equipped with an Air-Pouch at the tip. The catheter has one lumen that connects the Air-Pouch to the ICP-Monitor and one lumen that takes the guide wire for introduction.

## IAP-Catheter

The IAP-Catheter measures Intraabdominal Pressure with an Air-Pouch in the tip region of the catheter. It is compatible with all Spiegelberg ICP-Monitors. The IAP-Catheter is introduced like a naso-gastric tube.

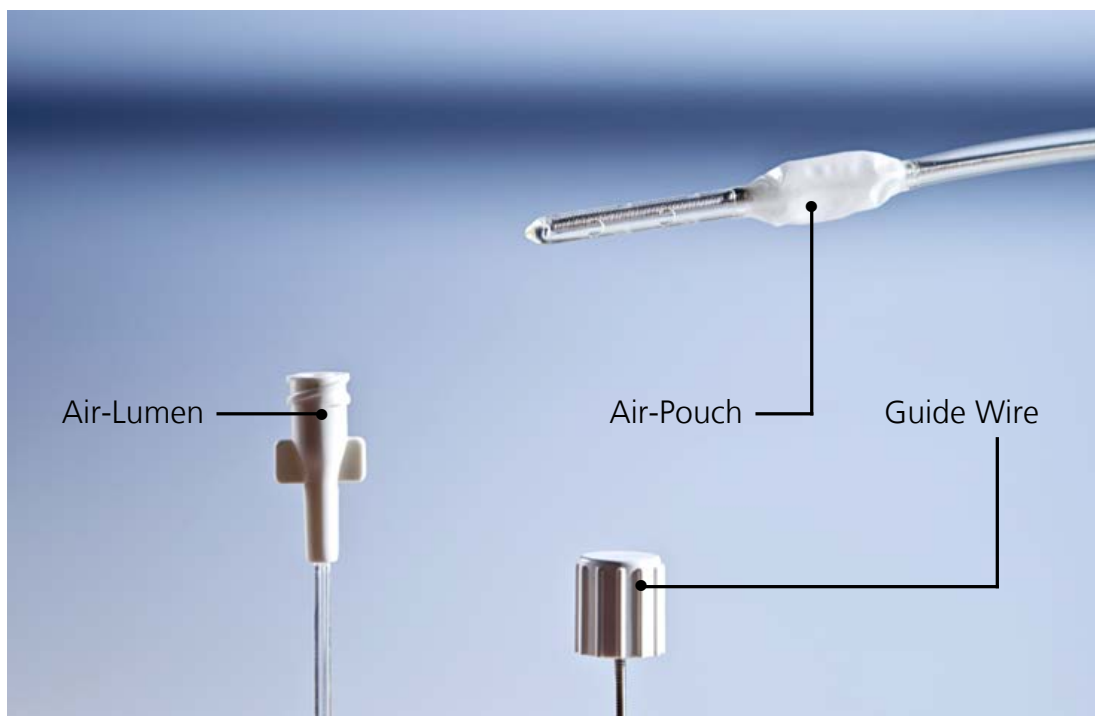

### Technical Information

Material	Polyurethane
Filling volume	0.05 - 0.1 cc
Outside Diameter	3.3 mm (10F)
Overall Length	2000 mm

### Ordering Information

<b>Article</b>	<b>REF</b>
IAP-Catheter	SND32.1.11

Spiegelberg GmbH & Co KG  
Tempowerkring 4  
21079 Hamburg  
Germany

Tel: +49 40 790178 0  
Fax: +49 40 790178-10

E-Mail: [info@spiegelberg.de](mailto:info@spiegelberg.de)  
WEB: [www.spiegelberg.de](http://www.spiegelberg.de)